



PRESSALIT®

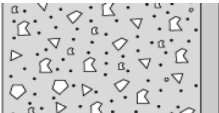
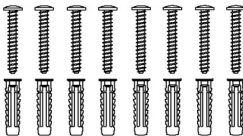

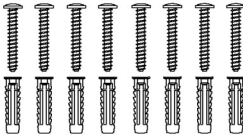
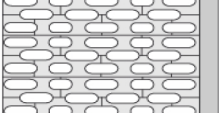
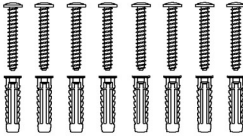
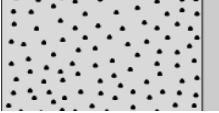
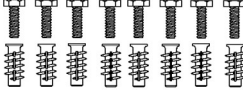
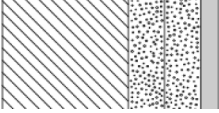
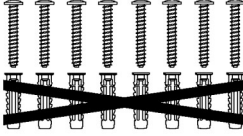
Product	PLUS handrail 1000 mm
Item number	RT104035
Cleaning	<p>Handrails ALU.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS:</p> <p>Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70%</p> <p>NON-DISINFECTANTS:</p> <p>Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS:</p> <p>Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>
Declaration of conformity	<p>CE UK MDR This product is in compliance with all relevant requirements of:</p> <p>- Regulation (EU) 2017/745 on medical devices.</p> <p>For the UK: - Medical Device Regulation (MDR) 2002 on medical devices.</p>



PRESSALIT®

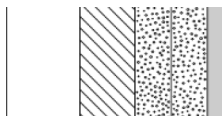
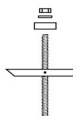

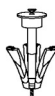
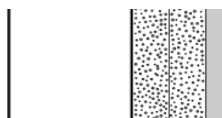
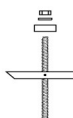

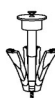
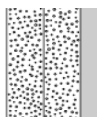

Product	PLUS handrail 1000 mm
Item number	RT104035

Installation in different wall types

	<p>Concrete >B25</p> <p>V8616 - Mounting kit, 8 screws Ø5 x 50 - plugs SX8, for brick wall. Must be ordered separately.</p>	
	<p>Solid Brick >= Mz12</p> <p>V8616 - Mounting kit, 8 screws Ø5 x 50 - plugs SX8, for brick wall. Must be ordered separately.</p>	
	<p>Solid Sand Lime Brick >= Hlz12</p> <p>V8616 - Mounting kit, 8 screws Ø5 x 50 - plugs SX8, for brick wall. Must be ordered separately.</p>	
	<p>Aerated Concrete P4</p> <p>V8666 - Mounting kit, 8 bolts M5 x 50 - plugs FTP K 6. To be ordered separately.</p>	
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8616 - Mounting kit, 8 screws Ø5 x 50 - plugs SX8. NB! Plugs are not needed. To be ordered separately.</p>	



PRESSALIT®

Product Item number	PLUS handrail 1000 mm RT104035			
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (8 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.			
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (2 mm) Wall fasteners (8 pcs) must be bought separately, for drywall it is possible to use e.g. fischer KD8 toggles or similar.			
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.			

pressalit.com