




PRESSALIT®

Product Item number	Mattress for MCT 2 and MCT 3 R842717378
Colour	VVS: 777542718 EAN: 5708590373881 HMI: 128440 Navy blue
Product description	 Mattress for MCT 2 and MCT 3 change trolleys, 695 x 1680 mm. Can also be used with MSCT 1 Mattress for MCT 2 and MCT 3 change trolleys, 695 x 1680 mm.
Materials	The mattress cover is manufactured from 100% polyester with polyurethane, approved in accordance with Oeko-Tex 100, Class I. It is fire retardant and impermeable to dust mites. The foam filling is flexible polyester-based polyurethane foam, which is also fire retardant and approved in accordance with Oeko-tex 100, Class I.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date.

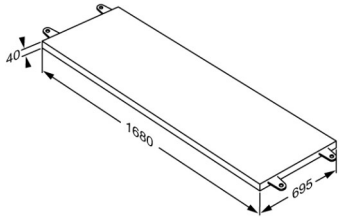

PRESSALIT®

Product	Mattress for MCT 2 and MCT 3
Item number	R842717378
Cleaning	<p>Mattress for CT and MCT. The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Sprit (husholdningssprit) Tvättsprit 70% Safe Zone Plus Selgiene Extreme Accel PREvention Wipes Lingettes</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>

pressalit.com

PRESSALIT®